

# CleanStream™ Dispenser System



**For Use with 1 or 2 dental units.**

## Directions for Use:

1. Add tap water to line marked "2 L".  
Screw on cap of bottle.
2. Cut off dotted corner of one CleanStream Evacuation System Cleaner Pouch. Pour □ contents of pouch completely into □ CleanStream Dispenser System.
3. Shake contents.
4. Clip atomizing adapters into their nest located in cap of bottle.
5. Attach saliva ejector to smallest of atomizing adapters.
6. Attach 1 or 2 high volume ejectors to respective adapter(s).
7. With vacuum pump on and handpiece valves open, aspirate CleanStream solution.
8. For first operatory cleansing, bottle must be in upright position.
9. For second cleansing lay bottle on its flat side.
10. After each cleansing, disconnect the handpiece.

## To reorder parts or entire system, refer to this Diagram:

### Tips

PN ED 598



PN 57643 SMALL



PN ED 597




PN 57644



Cap	PN 57645
O-Rings	PN 57647
Bottle	PN 57642

## Directions for multipurpose tip usage:

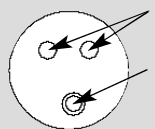

Clip tips directly into identical nests.

P/N ED 597



Use as adapter for post w/ Saliva Ejector

\*Note: Use w/ 2 HVE's

PN 57644



PN ED598



High Volume Evacuator



PN ED 597



Once attached these pieces stay together.

Adapter shown with Cap Overview on left.

PN 57643 Small



Saliva Ejector

