

Royal Dental Chair


Quick Reference Operator Control

Last-Position Operation (LP)

1. Push and Release LP. Chair moves from the current position to a pre-programmed "interim" position. Pushing the LP again will return the chair to the last position.

Programming LP "Interim"

1. Run chair to desired interim position.
2. Push and hold AP. While holding AP push and also hold LP until you hear 2 beeps (approx. 5 sec.). LP interim position is now set.

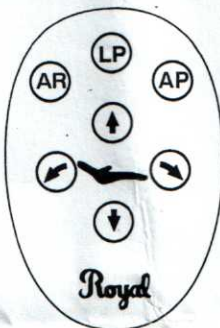

Auto-Position Operation (AP)

1. Push and release AP. Then push and release any manual control switch.
- NOTE: Each manual control switch can be programmed for a different position.

Programming Auto-Position

1. Run chair to desired position using manual controls.
2. Push and hold AP. While holding AP push and also hold preferred manual control switch until you hear 2 beeps (approx. 5 sec.). One Auto-Position (AP) is now set.

Auto-Return Operation (AR)


Push and release AR.
Chair goes to exit position.
(Not programmable)

Manual Controls

Base up/down.
Back up/down.


Optional 4 User Display

To change operator, push and hold AR. While holding AR, also push and release AP. Repeated pushes of AP will count through 4 different operators.

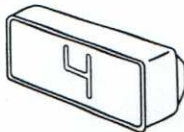

Chair Backrest Operation

(On select models.)

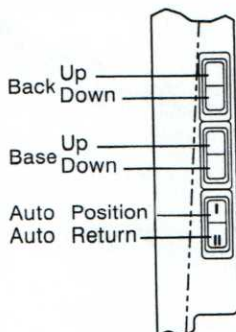

Royal Dental Manufacturing, Inc.
12414 Highway 99 Everett, Wa 98204-5548
Ph (425) 743-0988 Fax (425) 743-3588