


## Weekly Maintenance



### • Clean External Surfaces & Gaskets

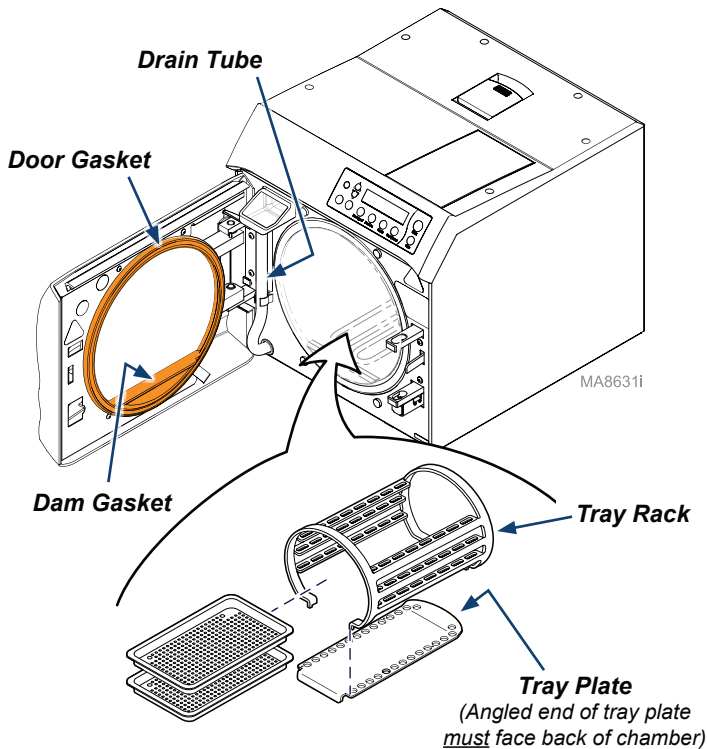

#### Caution

To prevent burns, allow unit to cool before cleaning gaskets and internal surfaces.

- Wipe external surfaces & gaskets with a soft damp cloth and mild soap / detergent.
- Examine door gaskets for damage.

### • Clean Internal Surfaces




#### Equipment Alert

Failure to change water may result in sterilizer malfunction.

- Drain water from the reservoir using drain tube on front of unit.
- Wash trays, rack, plate, and inside of chamber with mild soap or Speed-Clean and distilled water.
- Refill reservoir with clean, distilled water.

## Monthly Maintenance

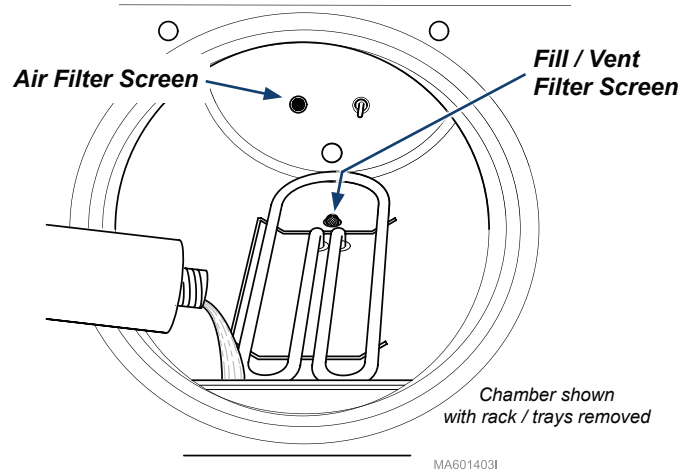

#### WARNING




Do not process instruments while flushing system.



#### Equipment Alert

Failure to flush system with **SpeedClean**, or use of **other** sterilizer cleaner(s), may result in the premature failure of sterilizer components.

### • Flush System

- With a cooled chamber, drain reservoir and fill with clean, distilled water. Add one ounce of SpeedClean sterilizer cleaner directly to the bottom of chamber.
- Run one **Pouches** cycle. 
- Press **Stop** button when Dry Cycle begins. 
- Drain reservoir and refill with clean, distilled water.
- Rinse by running one **Unwrapped** cycle. 
- Drain and refill reservoir with clean distilled water, then allow sterilizer to cool.

### • Clean Filter Screens (fill / vent , and air)



#### Equipment Alert

Do not operate sterilizer without filters in place.

- Allow sterilizer to cool, then remove trays / rack / plate.
- Refer to the illustration for location of filter screens.
- Pull / twist filter to remove. (Use pliers if necessary).
- Clean filters with SpeedClean and distilled water. (Replace filters if debris cannot be removed by cleaning).
- Rinse filters with distilled water.
- Install filters. (Press inward while twisting slightly).
- Install tray plate, rack, and trays.

### • Check Pressure Relief Valve\*

### • Remove / Clean / Inspect Gaskets\*


## To start a cycle...



### WARNING

If the Stop button is pressed before the Dry Cycle begins, the load is not sterile.

- Press the desired **Cycle Button**.  
(Refer to 'Cycle Parameters' for times / temperatures).
- Press the **Start** button.

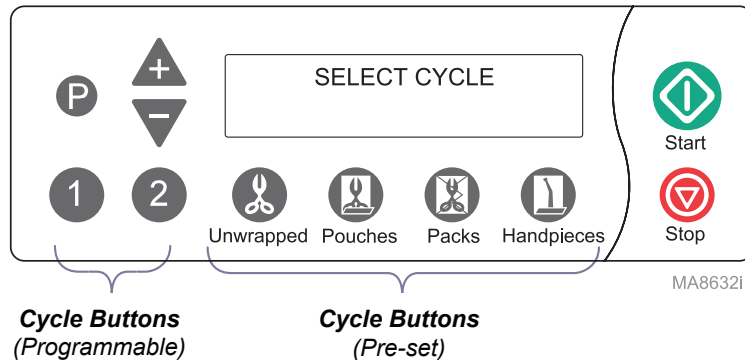

MA8632i

## To Program Cycle Buttons...



### Caution

Sterilization temperature can be adjusted from a minimum of 230°F (110°C) to a maximum of 275°F (135°C). Permitted temperature range for proper sterilization is 250° to 275°F (121° to 135°C). Temperatures set below 250°F (121°C) should not be used for sterilization, unless otherwise required by the device manufacturer. Temperatures below 250°F (121°C) are provided for disinfection only.

Step	Action	Description
1	Press desired button: 	This selects the button that will be programmed.
2	Press:	Program sterilization temperature (3 to 90 minutes)
	Adjust:	The "+" and "-" buttons adjust the temperature by 1° increments.
3	Press:	Program sterilization time
	Adjust:	The "+" and "-" buttons adjust the time by 1-minute increments.

Step	Action	Description
4	Press:	Program vent speed
	Adjust:	Pressing "+" changes setting: FAST Pressing "-" changes setting: SLOW
5	Press:	Program dry time (0 to 60 minutes)
	Adjust:	The "+" and "-" buttons adjust the time by 1-minute increments.
6	Press:	The display will show the new cycle parameters.

**NOTE:**  
Pressing the STOP button during this procedure will abort the changes, and revert to the original settings.  
The programmed settings are retained under Program # button <1> or <2>. Even if power is interrupted, or the unit is unplugged the setting will be retained.

### Cycle Parameters

	270°F (132°C) 27.1 psi (186 kPa) Sterilize: 3 min. Dry: 30 min.*	<ul style="list-style-type: none"> <li>Instruments loose on a tray.</li> <li>Open glass or metal canisters.</li> <li>Tubing not used in surgical procedures.</li> <li>Loose items manufacturers recommend for exposure at 270°F (132°C).</li> </ul> <p>Note: The sterility of unwrapped items is compromised on exposure to a non-sterile environment.</p>
	270°F (132°C) 27.1 psi (186 kPa) Sterilize: 5 min. Dry: 30 min.*	<ul style="list-style-type: none"> <li>Pouched or loosely wrapped instruments.</li> <li>Multiple layers of instruments separated by fabric.</li> <li>Wrapped trays of loose instruments.</li> <li>Tubing not used in surgical procedures.</li> <li>Wrapped items manufacturers recommend for exposure at 270°F (132°C)</li> </ul>
	250°F (121°C) 15 psi (104 kPa) Sterilize: 30 min. Dry: 30 min.*	<ul style="list-style-type: none"> <li>Textiles and surgical packs wrapped for sterilization.</li> <li>Items, except liquids, manufacturers recommend for exposure at 250°F (121°C) for 30 minutes.</li> </ul>
	270°F (132°C) 27.1 psi (186 kPa) Sterilize: 6 min. Dry: 30 min.*	<ul style="list-style-type: none"> <li>Dental handpieces (wrapped or unwrapped)</li> </ul> <p>Note: Verify acceptability of sterilization parameters with handpiece manufacturer.</p>
	230°F (110°C) to 275°F (135°C) 6 to 31 psi (41 to 214 kPa) Sterilize: 3 to 90 min. Dry: 0 to 60 min.	<ul style="list-style-type: none"> <li>Items appropriate for user's defined parameters.</li> </ul> <div> <b>Caution</b>            Temperatures below 250°F (121°C) should only be used for disinfection.         </div>

\* Dry time can be changed from 0 to 60 minutes. Refer to 'Cycle Operation' in the Installation / Operation Manual.