


**Figure 15**  
**Mounting Plate for**  
**Dual Stud Wall**  
**Mounts**

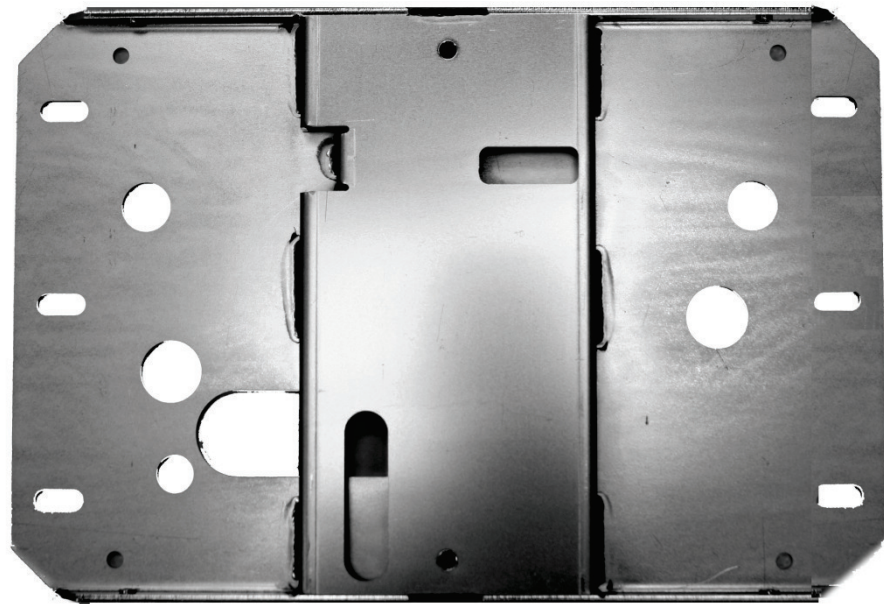

**Remove Control**  
**Unit Front Cover**

1. Open the shipping carton and locate the Control Unit in the first level of the carton.
2. Remove the Phillips screw from the front cover of the Control Unit.
3. Carefully remove the front cover.
4. Place the front cover and the screw in a safe location for later reassembly.

**Mount the**  
**Control Unit**

1. Select the 5/16-inch  $\frac{3}{4}$ -inch long socket cap screws and washers.
2. Put the upper screw and washer through the upper mounting hole of the Control Unit.
3. Place the Control Unit on the wall and loosely tighten the upper screw.
4. Put the lower screw and washer through the lower mounting hole of the Control Unit and loosely tighten. Be sure that the power wire extends through the opening at the bottom of the Control Unit.
5. Place a level on the Control Unit bearing parallel to the wall. Level the Control Unit.
6. Tighten the upper and lower screws.
7. After the Preva is installed, slide covers on the wall plate and screw on.