


A-dec 500 Wet/Dry Foot Control with Chip Blower/Accessory Button

# A-DEC 500 WET/DRY FOOT CONTROL WITH CHIP BLOWER/ACCESSORY BUTTON

- Modulates variable speed operation with foot pressure applied in any position
- Use for air turbine, electric handpiece, A-dec camera, as well as ultrasonic instruments
- Dark Surf or Gray 3
- For use with A-dec 500 delivery systems only

PRODUCT DESCRIPTION	PART NUMBER SURF	PART NUMBER GRAY 3	PRICE
A-dec 532/533 foot control 4\mathbb{A}" (1,321 mm)	38.1805.00	38.1829.00	⊠380.00
Tubing replacement kit	90.1104.00	90.1104.01	⊠40.00
A-dec 541/542 foot control, 16\(\tilde{\text{2}}\)(4,876 mm)	38.1805.01	38.1829.01	⊠380.00

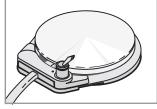

A-dec 300 Wet/Dry Foot Control

# A-DEC 300 WET/DRY FOOT CONTROL

- Variable speed operation with foot pressure applied in any position
- Toggle operates On/Off control for coolant water
- Standard tubing: 7⊠(2,134 mm) Gray 3
- Use for both the air turbine or electric handpiece

1	PRODUCT DESCRIPTION	PART NUMBER	PRICE
	A-dec 300 foot control Gray 3	38.1821.00	⊠300.00
	A-dec 300 euro conversion kit	38.0825.00	⊠59.00



Foot Control

# FOOT CONTROL III

- Variable speed operation with foot pressure applied in any position
- For use with Cascade or Radius units equipped with the Master Series touchpad option
- Standard tubing: 7⊠(2,134 mm) Dark Surf

PRODUCT DESCRIPTION	PART NUMBER	PRICE
Cascade foot control	38.1700.00	⊠280.00




Wet/Dry Foot Control

# WET/DRY FOOT CONTROL III

- Variable speed operation with foot pressure applied in any position
- Toggle operates On/Off control for coolant water
- Standard tubing: 7⊠(2,134 mm) Dark Surf
- Use for both the air turbine or electric handpiece

PRODUCT DESCRIPTION	PART NUMBER	PRICE
Dark Surf	38.1603.00	⊠360.00
Euro wet/rdy, Dark Surf	38.1603.01	⊠360.00
Foot control III euro conversion kit	38.1764.00	⊠59.00

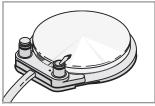

Wet/Dry Foot Control with Chip Blower/Accessory Button

# WET/DRY FOOT CONTROL III WITH CHIP BLOWER/ACCESSORY BUTTON

- Variable speed operation with foot pressure applied in any position
- Toggle operates On/Off control for coolant water
- Button operates chip blower valve or air-actuated electric switch for an ultrasonic instruments
- Standard tubing: 7⊠(2,134 mm) Dark Surf

PRODUCT DESCRIPTION	PART NUMBER	PRICE
Chip blower, Dark Surf	38.1601.00	⊠380.00
Euro wet/dry/chip blower, Dark Surf	38.1601.01	⊠380.00
Foot control III euro conversion kit	38.1764.00	⊠59.00

### FOOT CONTROL I VALVE SERVICE KIT

90.0010.00



• For use with Foot Control I Valve (38.0659.00)

# FOOT CONTROL III VALVE

38.0761.00

⊠130.00



- With relay
- Replaces 38.0607.01 and 38.0659.00
- Not for A-dec 500 foot control

# LEVER FOOT CONTROL SERVICE KIT

90.0313.00

⊠100.00



• For use with Lever Foot Control Valve (38.0201.00)

### FOOT CONTROL II VALVE SERVICE KIT

90.0312.00



• For use with Foot Control II Valve (38.0607.01)

# FOOT CONTROL III VALVE SERVICE KIT

90.0593.00



- For use with Foot Control III Valve (38.0761.00)
- Also use with A-dec 500 foot control

### **SERVICE TIP**

Water leaks from handpieces after removing foot from foot control

Some water leaks can be caused by pinched foot control tubing, which blocks the exhaust exit out of the foot control and traps the water coolant signal.

Reminder: Troubleshooting an Active Foot Control

Before removing the foot control disc, move the master On/Off toggle to the OFF position and bleed the system of air pressure. This prevents the foot control stem from being forcefully ejected from the Foot Control II valve.

If the air supply to the foot control must be active while testing, before moving the master On/Off toggle to the ON position, ensure that the stem assembly in the foot control valve cannot be ejected.

65


# WET/DRY TOGGLE (3-WAY) VALVE



### CHIP BLOWER BUTTON WITH 2-WAY VALVE

38.0610.00 Surf 6

**⊠75.00** 

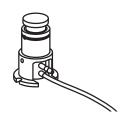

For use with Foot Controls


 I, II, and III as a chip blower
 button

# CHIP BLOWER/SCALER/ACCESSORY BUTTON WITH 3-WAY VALVE



- Use with A-dec 500 foot controls as chip blower/ Accessory button
- Use with Foot Controls I, II and III as a chip blower button

# **MICRO VALVE**

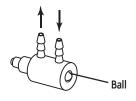

- 2-way (brass ball) for chip air
- 3-way (stainless steel ball) for wet/dry or scaler

# FOOT CONTROL I AND II ONE-HOLE (RIGHT) HOUSING

38.0320.02 Surf 6

⊠17.25



 Use with chip blower, ultrasonic instrument, wet/ dry or foot-activated water dispensing system

# FOOT CONTROL III ONE-HOLE (LEFT) HOUSING

38.0763.00 Surf 6

 $\boxtimes 17.25$ 



Use with Foot Control III with wet/dry toggle

# FOOT CONTROL HOUSING, TWO-HOLE



 Use with wet/dry with chip blower, wet/dry with ultrasonic instrument, or ultrasonic instrument with chip blower

USA and Canada Phone: 1.800.547.1883 : 24-hour Fax: 1.503.538.0276 | International Phone: 1.503.

International Phone: 1.503.538.7478 : 24-hour Fax: 1.503.538.5911

# Components

# **TOGGLE KIT**

38.0075.03 Surf 6 ⊠12.00 38.0075.04 Gray 3 ⊠12.00



• Includes pin, toggle, and retaining spring

# **RUBBER BUMPER PAD (FEET)**

⊠15.00 017.007.01 (25)



• 1/2" diameter bumper with self-adhesive back

# FOOT CONTROL I RETAINING RING

22.0120.00 Gray ⊠24.00



Includes 4 screws 6 x 1/2"

# **FOOT CONTROL COVER**

22.0110.00 ⊠45.00




- Fits all A-dec foot controls
- · Chrome plated

### **FOOT CONTROL TUBING** LENGTH NO SHEATH **GRAY 3** PRICE DARK SURF Tubing Without Sheath 10' (3m) 024.151.01 ⊠30.00 50' (15m) 024.151.03 ⊠125.00 Foot Control III Tubing | 10' (3m) 024.179.01 ⊠40.00 50' (15m) 024.179.03 ⊠160.00 A-dec 500 Foot Control Tubing 10' (3m) 024.202.01 024.229.01 ⊠40.00 50' (15m) 024.202.02 024.229.02 ⊠160.00 A-dec 300 Foot Control Tubing 10' (3m) 024.214.01 ⊠40.00 024.214.02 50' (15m) ⊠160.00

# FOOT CONTROL II AND III RETAINING RINGS

38.0237.00 Black

⊠6.50

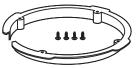

- For Foot controls II and III
- Includes 4 screws 6 x 1/2"

USA and Canada Phone: 1.800.547.1883 : 24-hour Fax: 1.503.538.0276 International Phone: 1.503.538.7478 : 24-hour Fax: 1.503.538.5911