



## 1. EzRay Control Panel Extension kit (Optional Purchase)

### Contents

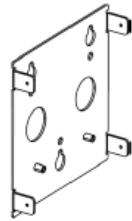

BRACKET BOTTOM

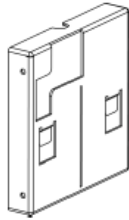


BRACKET TOP



PLASTIC COVER



EXTENSION CABLE  
(15m)

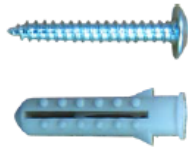

SCREW ANCHOR SET x 3  
(M4 x 30)



TRUSS BOLT x 6  
(M4 X 8)




CABLE CLAMP x 2




INSTALLATION  
GUIDE

### Overview

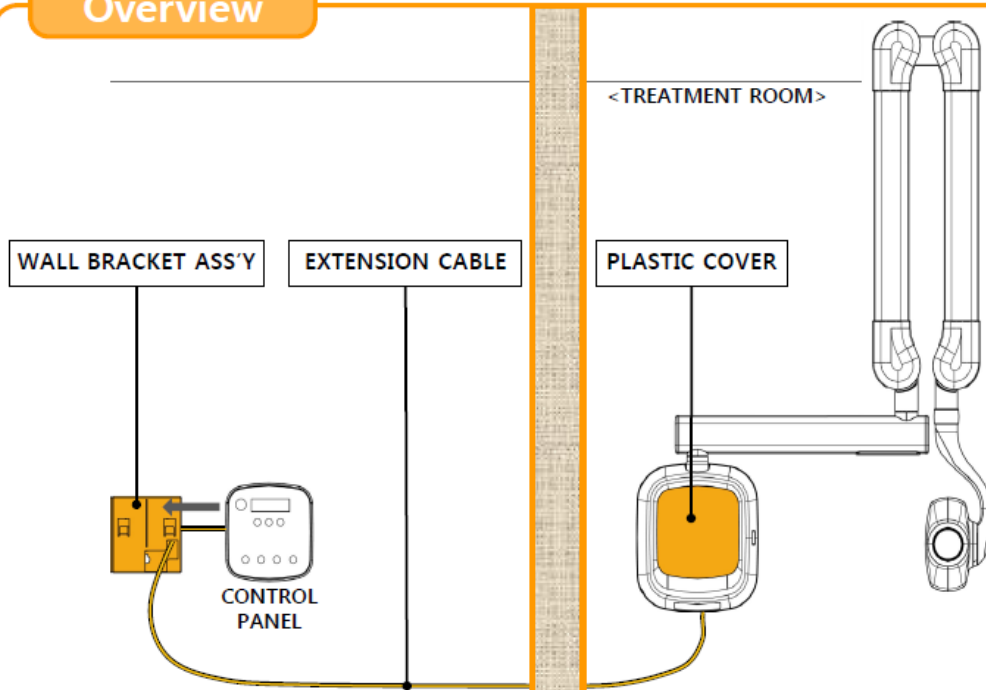

## 2. Mounting Options (NOTE: A custom backing is needed for a “Pass-thru Cabinet” installation, Single Stud Mount is currently under development)

### Mounting Type

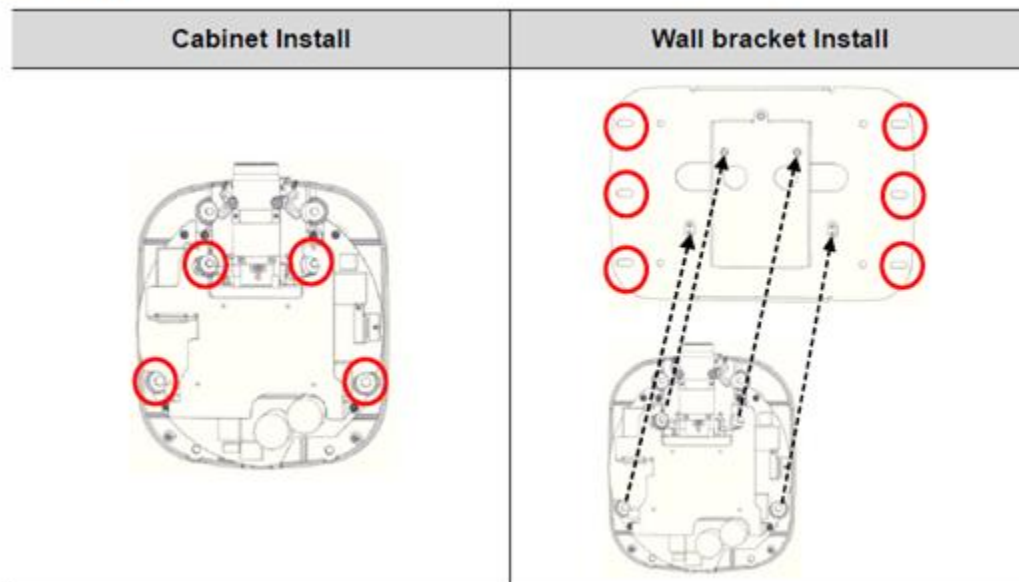

## 3. How to shut-off exposure switch on the control panel

- Press circled two buttons simultaneously (Enters admin mode).
- Use up, down buttons to go to setting #6, change control panel expos setting to off (press upper circled button to change).
- Press bottom right circled button for 5 seconds to “save” the setting.

