
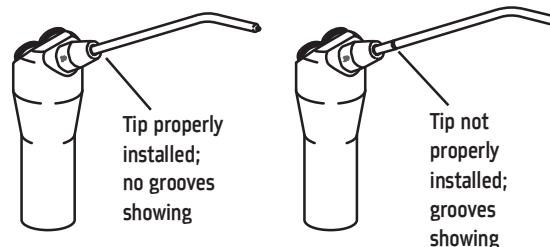


## USE ONLY A-DEC SYRINGE TIPS

Use of syringe tips manufactured by a company other than A-dec may result in damage to the O-rings inside the syringe tip retainer assembly. Damage to these O-rings can result in syringe tip ejection. Many disposable syringe tips have sharp corners on the water tube, which cut the O-ring sealing the water from the air in the tip.

To help you easily identify genuine A-dec syringe tips, a distinguishing "a" mark is located above three grooves on each tip. A-dec syringe tips also have a bright nickel chrome finish.

The nickel chrome finish on A-dec's syringe tips is not affected by cleaning solutions containing ammonia or amines. While A-dec tips are not affected by these solutions, we discourage the use of these chemicals when cleaning A-dec or any other manufacturers' syringe tips since tips from different manufacturers are often mixed together in the dental office.

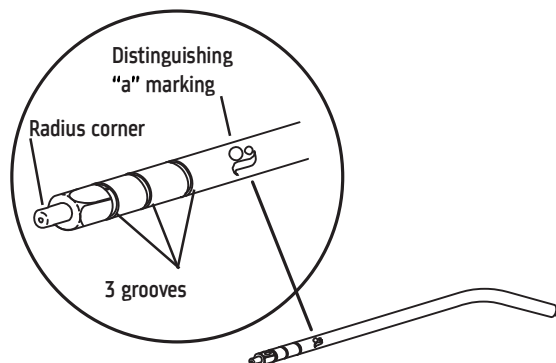

### STANDARD SYRINGE TIP

23.0872.01 (5) \$41.50


23.0872.02 (25) \$175.00



- Autoclavable
- Standard size 5<sup>1</sup>/<sub>4</sub>" (83 mm)
- Specially designed for A-dec syringes

### SYRINGE TIP RETAINER NON-LOCKING KIT

23.1112.00 Non-locking kit \$19.75

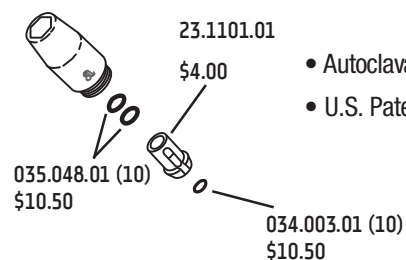


- Autoclavable
- U.S. Patent 5,433,485

### SYRINGE TIPS

23.0464.00 Long (1) \$30.00

23.0465.00 Short bend (1) \$15.00

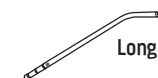

Long




Short bend

- Autoclavable
- Long size 5<sup>1</sup>/<sub>4</sub>" (133 mm)
- Short bend 3<sup>1</sup>/<sub>4</sub>" (83 mm)
- Specially designed for A-dec syringes


### SYRINGE TIP RETAINER LOCKING KIT

23.1090.00 Locking kit \$35.00

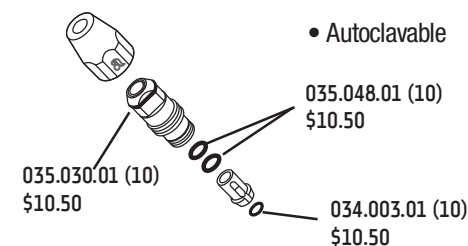


- Autoclavable

## AUTOCLAVABLE TRADITIONAL SYRINGE (beginning October 1995)

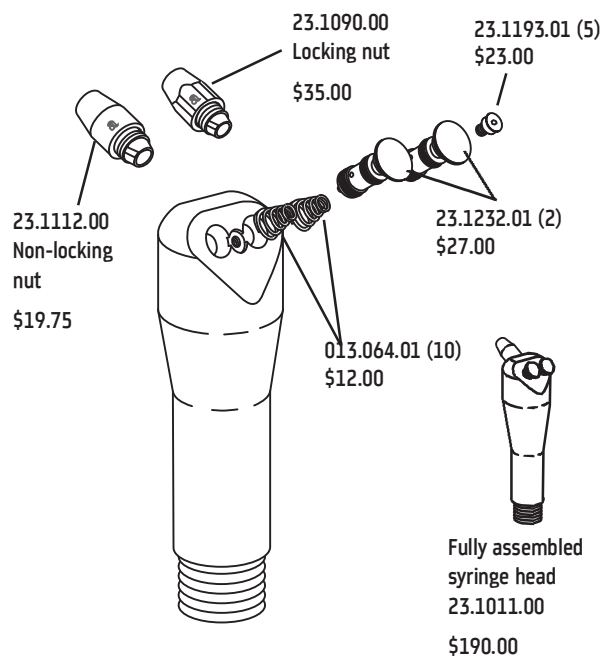

## AUTOCLAVABLE TRADITIONAL SYRINGE (prior to October 1995)


Service Parts Only

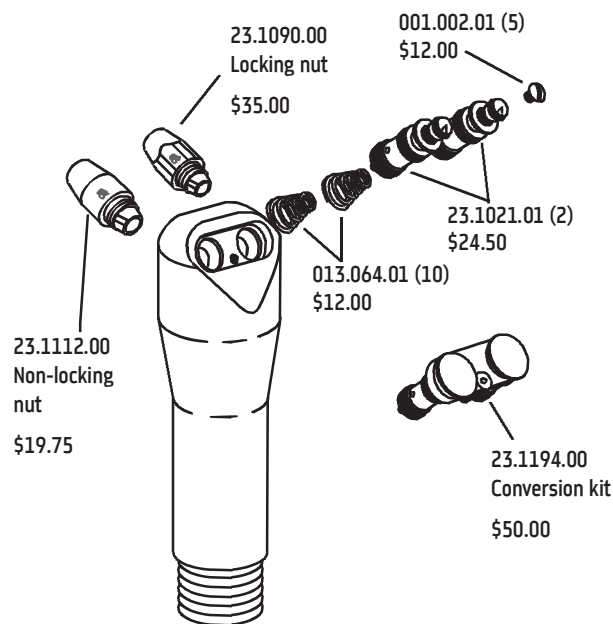

## TERMINAL BARB ASSEMBLY NON-QUICK-DISCONNECT

23.1110.00 2-barb


☒\$25.00

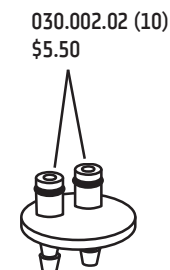


## TERMINAL BARB ASSEMBLIES QUICK-DISCONNECT

23.1067.00 2-barb


☒\$50.00

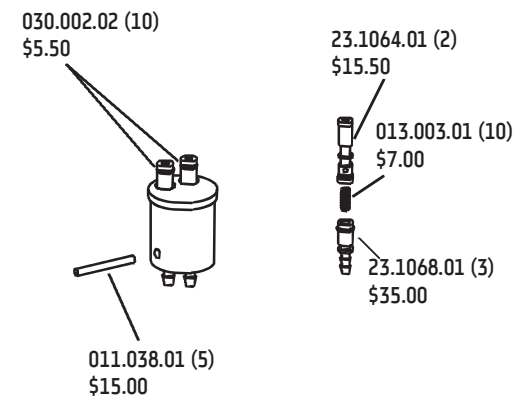

## STANDARD SYRINGE KITS


	SYRINGE ASSEMBLIES (Include syringe, tubing and flow control clamp)	PART NUMBER SURF	PART NUMBER GRAY 3	PRICE
	With Terminal Barb Assembly 7/8(2,134 mm)	23.1150.00	23.1302.00	☒\$260.00
	With Terminal Barb Assembly 12/32(3,657.6 mm)	23.1150.12	23.1311.00	☒\$260.00
	With QD Terminal Barb Assembly 7/8(2,134 mm)	23.1195.00	23.1303.00	☒\$290.00

## AUTOCLAVABLE CONTINENTAL SYRINGE (Beginning October 1995)

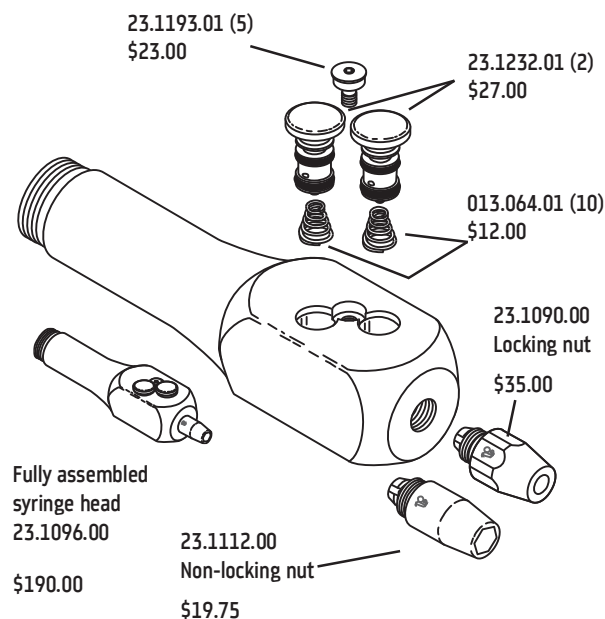

## AUTOCLAVABLE CONTINENTAL SYRINGE (prior to October 1995)


Service Parts Only

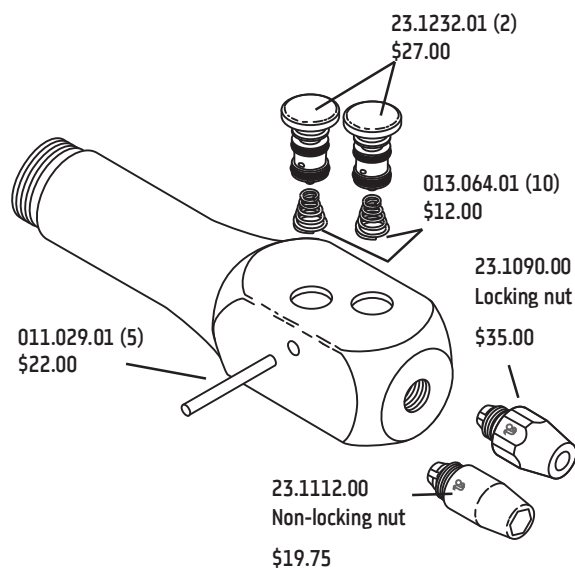

## TERMINAL BARB ASSEMBLY NON-QUICK-DISCONNECT

23.1110.00 2-barb

☒25.00

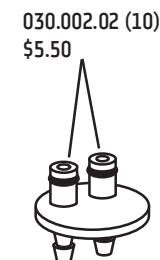


## TERMINAL BARB ASSEMBLIES QUICK-DISCONNECT

23.1067.00 2-barb


☒50.00

23.1066.00 3-barb, circulating

☒125.00

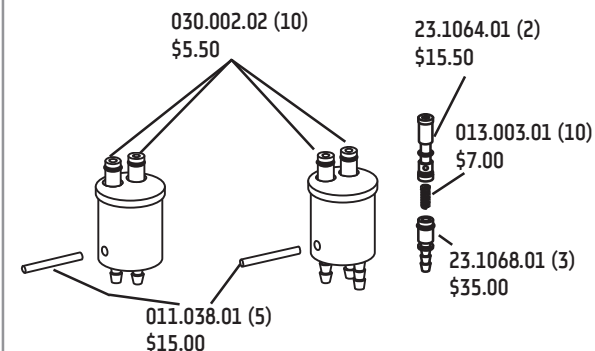

## AUTOCLAVABLE CONTINENTAL SYRINGE

	SYRINGE ASSEMBLIES (Include syringe and tubing)	PART NUMBER SURF	PART NUMBER GRAY 3	PRICE
	With Terminal Barb Assembly 84" (2,134 mm)	23.1170.00	23.1304.00	☒260.00
	With QD Terminal Barb Assembly 84" (2,134 mm)	23.1247.00	23.1305.00	☒290.00

**BUTTON VALVES**

23.1232.01 Buttons (2)

☒27.00




- Requires spring 013.064.01 (10) ☒12.00

**AUTOCLAVABLE SYRINGE VALVE KIT**

23.1234.00

☒33.00

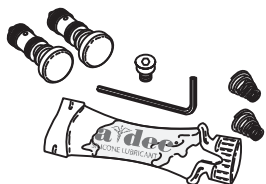

- For syringes beginning October 1995
- Includes button valves, springs, retaining screw, hex key, and silicone lubricant

**BUTTON VALVES**

23.1021.01 Buttons (2)

☒24.50




- Requires spring 013.064.01 (10) ☒12.00

**AUTOCLAVABLE SYRINGE TIP CONVERSION KIT**

23.0012.00

☒67.00

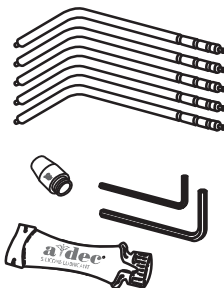


- Converts a locking syringe tip into a non-locking

**2-BUTTON CONVERSION KIT**

23.1194.00

☒50.00

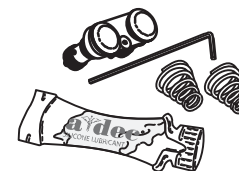

- Converts from single autoclavable button to two autoclavable buttons
- Includes: (2) valve assemblies, (1) spacer, screw, (2) springs, silicone lubricant, and hex key

**AUTOCLAVABLE QD SYRINGE SERVICE KIT**

23.1235.01

☒93.00



- Includes syringe tip, screw, pin, syringe nut with spacer, springs, button valves, stem assemblies, terminal barbs, silicone lubricant and O-rings

**SERVICE TIPS**

Use only A-dec silicone lubricant 98.0090.01


Using petroleum-based lubricants (jellies, creams, or lotions) causes the O-rings in the syringe to swell and/or deteriorate, resulting in syringe malfunction or failure.

Use only A-dec syringe tips

Using disposable syringe tips or tips manufactured by a company other than A-dec may result in damage to the O-rings inside the tip retainer nut assembly. Damaged O-rings may result in syringe tip ejection.


Identify genuine A-dec tips by the ☒☒ located above three grooves on each tip.

## EZ-CLEAN SYRINGE SERVICE PARTS ONLY

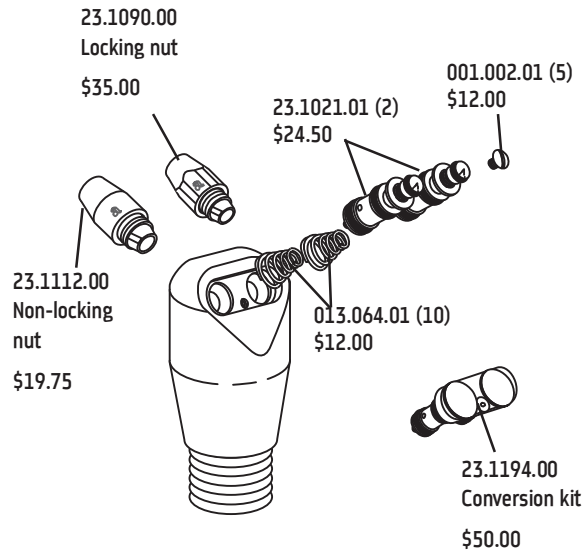

## CONTINENTAL PENCIL SYRINGE


Service Parts Only

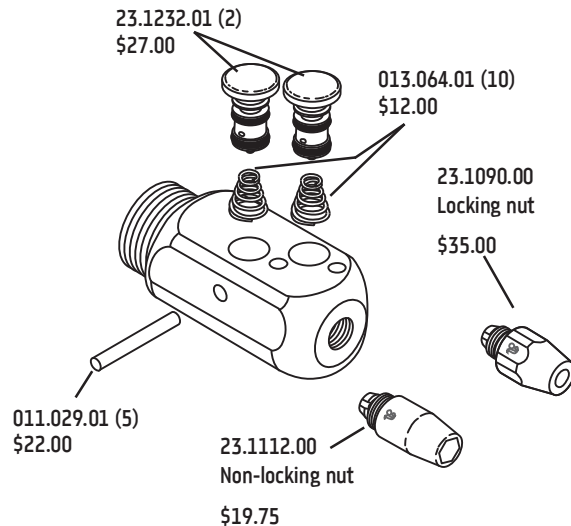

## TERMINAL BARB ASSEMBLIES

23.0760.00 Non-circulating


☒70.00

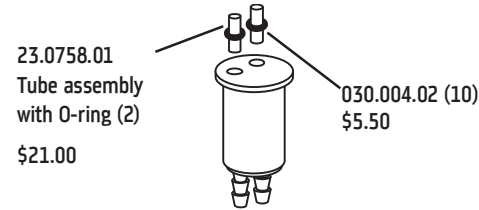

## EZ-CLEAN SYRINGE CONVERSION KIT

23.1060.00

☒250.00

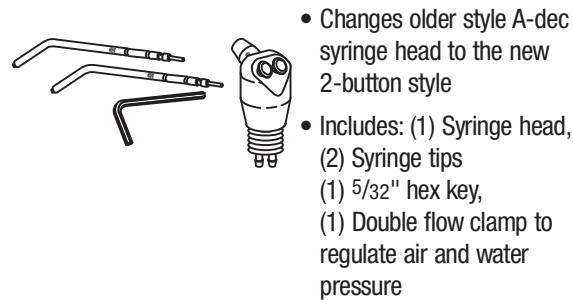

## BUTTON VALVE FOR EZ-CLEAN SYRINGE

23.1021.01 (2)

☒24.50




## SERVICE TIP

A burst of water from the tip when the syringe air button is pressed could be the result of a worn or damaged O-ring (034.003.01) in the tip nut spacer. This O-ring prevents water-air crossover and should be replaced once a year.


## LEVER SYRINGE

Service Parts Only

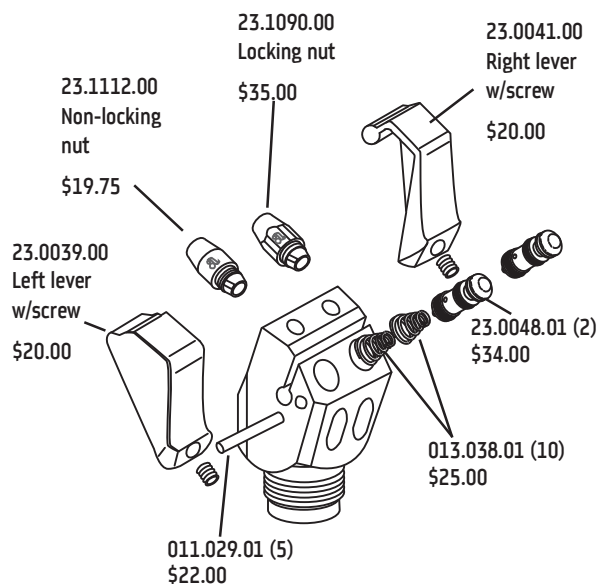

## SOFT-TOUCH SYRINGE

Service Parts Only

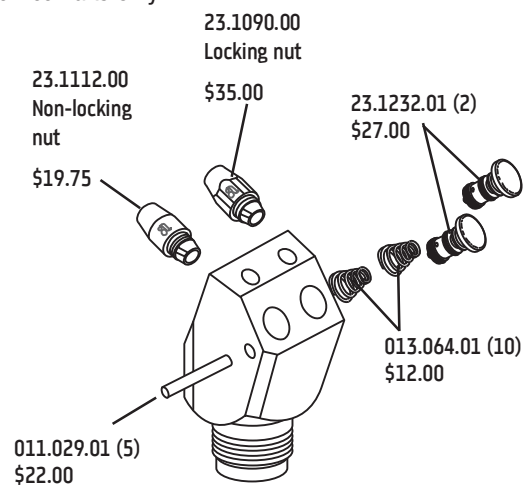

## BUTTON VALVE FOR LEVER SYRINGE

23.0048.01 (2)

☒34.00



- Requires spring 013.064.01 (10) ☒12.00

## LEVER SYRINGE VALVE KIT

23.1093.00

☒40.00



- To replace valves on lever syringes
- Includes button valves, springs, and silicone lubricant

## AUTOCLAVABLE SYRINGE WITH QD SERVICE KIT

23.1235.01

☒93.00



- Includes syringe tip, screw, pin, syringe nut with spacer, springs, stem assemblies, valve assemblies terminal barbs, silicone lubricant, and O-rings


## LEVERS

23.0039.00

☒20.00

23.0041.00


☒20.00

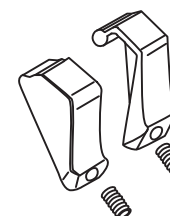

- Includes adjustment screw


## UNI-FLO SYRINGE

Service Parts Only

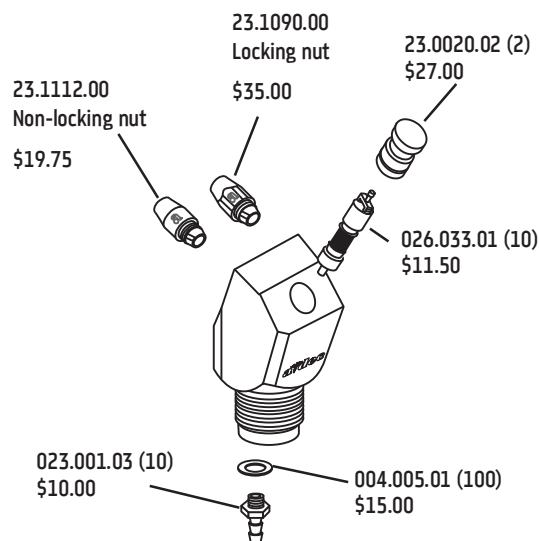

## TRI-FLO SYRINGE

Service Parts Only

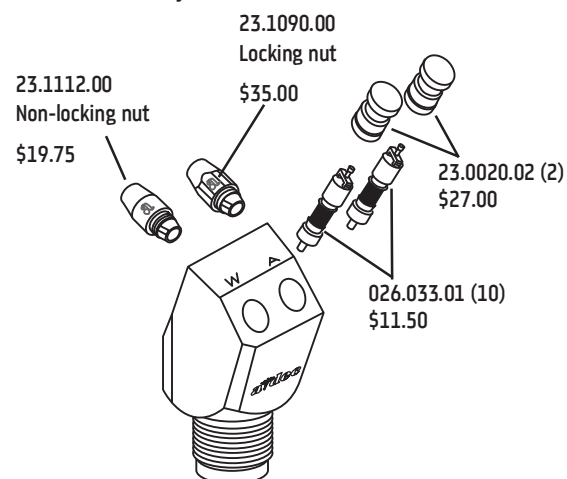

## BUTTON VALVE

23.0020.02 (2)

☒27.00



## SYRINGE BUTTON DIAPHRAGMS

23.0220.02 (10)

☒19.00




- For use on 23.0020.02 syringe button assy

## UNI-FLO AND TRI-FLO SYRINGE VALVE CORE SERVICE KIT

23.0018.00

☒21.50

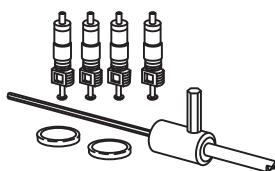


- Includes four valve cores, two syringe button diaphragms, and valve core removal tool

## UNI-FLO AND TRI-FLO SYRINGE SERVICE KIT

23.0019.00

☒47.00

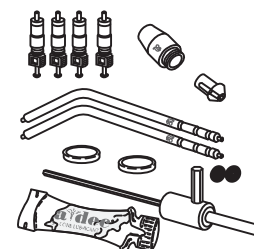


- Includes two syringe tips, nut, four valves, two diaphragms, two O-rings, valve core removal tool, and lubricant

## UNI-FLO SYRINGE SERVICE KIT

23.0027.00

☒35.00

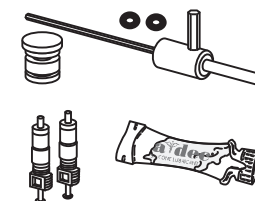


- Includes one button assembly, two valves, two O-rings, valve core removal tool, and lubricant

## A-DEC SYRINGE TOOL

98.0010.00

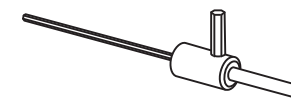

☒14.25



A-dec's syringe tool incorporates:

- A valve core removal tool
- A 3/32" hex wrench for syringe flow adjustment
- A 5/32" hex wrench for removing and installing the syringe tip insert

## SYRINGE TUBING

		PART NUMBER SURF	PART NUMBER GRAY 3	PRICE
	Straight 3-hole 7/8(2,134 mm)	024.155.02	024.223.02	¥30.00
	Straight 3-hole 1 1/2(3,657.6 mm)	024.155.12	024.223.12	¥30.00
	Assembly 7/8(2,134 mm) (Includes collet, bushing, retainer, retainer spring, and uni-clamps)	23.1208.00	23.1208.01	¥39.50

## REPLACEMENT HANDLE

23.1015.00 ¥15.00



- Use for:  
Autoclavable, and  
Continental autoclavable  
syringes

## REPLACEMENT HANDLE

23.0757.00 Handle ¥60.00

23.0768.00 Bushing ¥23.00



Collet

- Use for :  
Continental pencil, EZ  
Clean, Lever, Soft-Touch,  
Uni-Flo and Tri-Flo
- Use collet with  
Uni-Flo, syringe

## SYRINGE STRAIN RELIEF REPAIR KIT

90.0563.00 Repair kit ¥12.00

010.073.01 Spring (5) ¥18.00



Springs



O-ring

- Includes (1) O-ring and  
(2) retainer springs



## A-DEC WARM WATER SYRINGE ASSEMBLIES




- Include warm water syringe and tubing assembly, warm water circuit board, and attaching hardware



### WARM WATER SYRINGE

PRODUCT DESCRIPTION	USAGE	☒795.00		☒365.00		☒180.00
		Full Assy		Tubing		Syringe Head
		SURF	GRAY 3	SURF	GRAY 3	
SYR,WM,TRD,84"/56"SHTH	532, 541, 542, 332, 334	23.1272.01	23.1290.01	23.1284.00	23.1301.01	23.1278.00
SYR,WM,CON,84"/56"SHTH	533, 333, 335	23.1273.01	23.1291.01	23.1284.00	23.1301.01	12.1277.00
SYR,WM,TRD,120"/56"SHTH	551, 351	23.1274.01	23.1306.01	23.1279.01	23.1307.01	12.1278.00
SYR,WM,TRD,84"/62"SHTH	545	23.1283.01	23.1294.01	23.1285.00	23.1299.01	12.1278.00

### LOW VOLTAGE WATER HEATER SERVICE PARTS

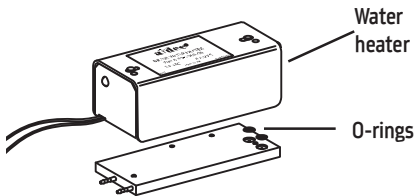

Tubing assembly 5' Dark Surf

DESCRIPTION	PART NUMBER	PRICE
Water heater, low voltage	40.1060.00	☒550.00
O-ring (10)	033.003.01	☒27.00
Heated tubing assembly with terminal (shown above)	23.1198.00	☒99.00
Heated tubing assembly, no terminal	024.170.01	☒73.00