

Classic Series™ & PowerAir™ Compressor Pre-Installaion Guide

Applies to Models: CL21, CL22, CL32, CL52 P21, P22, P32, P52, P72

Classifications: Class 1, Type B Applied Part

Classic Series™ Specifications

Model	# of Users	CFM @ 80 PSI	Total HP	Tank Capacity	Voltage	Total Amps	Breaker Size (Amps)	Sound Level (dbA)	Dimensions H x W x D (IN.)	Product Weight (lbs.)
CL21	1 - 3	5.3	1.0	10	115	12.4	20	68	29 x 26 x 24	202
CL22	1 - 3	5.3	1.0	10	208-230	6.0	20	68	29 x 26 x 24	202
CL32	3 - 5	10.6	2.0	20	208-230	12.0	20	71	31 x 33 x 26	329
CL52	5 - 7	15.9	3.0	32	208-230	18.0	30	75	32 x 41 x 27	460

PowerAir™ Specifications

Model	# Users	CFM @ 80 PSI	Total HP	Tank Capacity	Voltage	Total Amps	Breaker Size	Sound Level (dBA)	Sound Level (dBA) With Quiet Cover	Dimensions H x W x D (IN.)	Dimensions With Sound Cover H x W x D (IN.)	Product Weight (LBS.)
P21	1 - 3	5.2	1.5	10	115	16.0	20	65	60	27 x 26 x 24	32 x 26 x 25	225
P22	1 - 3	5.2	1.5	10	208-230	7.0	20	65	60	27 x 26 x 24	32 x 26 x 25	225
P32	3 - 5	7.8	2.25	20	208-230	11.0	20	67	62	29 x 33 x 26	33 x 33 x 28	295
P52	5 - 7	10.6	3.2	20	208-230	20.0	30	70	65	31 x 33 x 26	33 x 33 x 28	345
P72	7 - 10	15.9	4.8	32	208-230	29.0	40	73	68	33 x 40 x 27	39 x 40 x 30	425

Compressor Site Requirements

Electrical Models

Supply	CL21	CL22	CL32	CL52	P21	P22	P32	P52	P72
Voltage	115	208-230	208-230	208-230	115	208-230	208-230	208-230	208-230
Single Phase					60 Hz				
Disconnect Switch Box									
- Fused									
- User Supplied	20 Amps	20 Amps	20 Amps	30 Amps	20 Amps	20 Amps	20 Amps	30 Amps	40 Amps
- Box(es) must be located within	-			-				-	
3 ft. of compressor									

Plumbing

Fresh Air Intake Line						
Туре	PVC					
Line Sizing	1" Minimum (Provided by Contractor)					
Compressor Connections	1" MPT or 1" PVC Slip Connectors					
Drain						
Туре	Poly Tube (Provided with Install Kit)					
Delivery Line						
	3/8" to 3/4" Copper, as site layout dictates. Terminate in mechanical room with 1/2" FPT. Flex hose provided for service					
Main Trunk Line	connection. Trunk line to slope 1/4" / 10ft toward source compressor.					
Operatory Branch Lines	3/8" to 3/4" Copper - Same size as trunk line. Terminate 1/2" or 5/8" FPT at Operatory junction box location.					

Environmental

Equipment Room Ambient	40° - 104° Fahrenheit / 4° - 40° Celsius
Temperature - Operational	40° - 104° Fahrenheit / 4° - 40° Celsius

Installation Kit Provided with Compressors



Exhaust Silencer

Indicator for Coalescing Filter



Fresh Air Intake Hoses Length 10'











Feet with Hardware

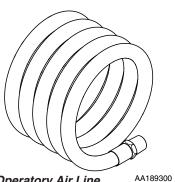


Operatory Air Line Length - 6'

Compressor Site Layout

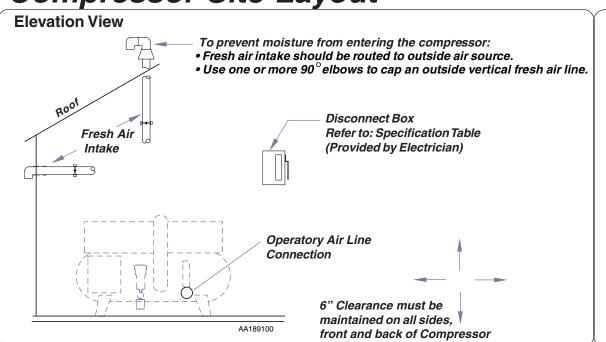

OR

• Fresh air intake can be set-up to pull air from a drop celing or another room.

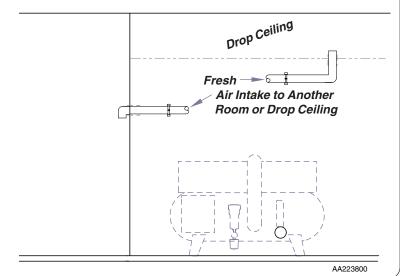

Sample Plumbing Layout

AIR LINE SIZING CHART						
Compressor CFM	Main Header Line					
@ 80 psig	Size					
5 - 10	3/8" ODT					
10 - 20	1/2" ODT					
20 - 30	5/8" ODT					
30 - 40	3/4" ODT					

